GM130 Rabbit mAb [7lqZ]

Cat NO. :A71733

Information:

Applications	Reactivity:	UniProt ID:	MW(kDa)	Host	Isotype	Size
WB IHC ICC/IF IP	Human,Mouse,R	Q08379	130kDa	Rabbit	lgG	50ul,100ul,200ul
	at Monkey Cow					
	Dog					

Applications detail:

Application	Dilution			
WB	1:1000-2000			
IHC	1:100			
ICC/IF	1:100			
The optimal dilutions should be determined by the end user				

Conjugate:

UnConjugate

Form:

Liquid

sensitivity

Endogenous

Purification:

Affinity-chromatography

Specificity:

Antibody is produced by immunizing animals with A synthesized peptide derived from human GM130

Storage buffer and conditions:

Antibody store in 10 mM PBS, 0.5mg/ml BSA, 50% glycerol (buffer) .

Shipped at 4°C. Store at-20°C or -80°C.

Products are valid for one natural year of receipt. Avoid repeated freeze / thaw cycles.

Tissue specificity:

Subcellular location:

Golgi apparatus, cis-Golgi network membrane,Peripheral membrane protein,Cytoplasmic side. Endoplasmic reticulum-Golgi intermediate compartment membrane,Peripheral membrane protein,Cytoplasmic side.

Introduction: WB: Western Blot IP: Immunoprecipitation IHC: Immunohistochemistry ChIP: Chromatin Immunoprecipitation ICC/IF: Immunocytochemistry/ Immunofluorescence F: Flow Cytometry


Cross Reactivity: H: human M: mouse R: rat Hm: hamster Mk: monkey Vir: virus MI: mink C: chicken Dm D. melanogaster X: Xenopus Z: zebrafish B: bovine Dg: dog Pg: pig Hr: horse

For Research Use Only. Not For Use In Diagnostic Procedures.

Peripheral membrane component of the cis-Golgi stack that acts as a membrane skeleton that maintains the structure of the Golgi apparatus, and as a vesicle thether that facilitates vesicle fusion to the Golgi membrane (Probable) (PubMed:16489344). Required for normal protein transport from the endoplasmic reticulum to the Golgi apparatus and the cell membrane (By similarity). Together with p115/USO1 and STX5, involved in vesicle tethering and fusion at the cis-Golgi membrane to maintain the stacked and inter-connected structure of the Golgi apparatus. Plays a central role in mitotic Golgi disassembly: phosphorylation at Ser-37 by CDK1 at the onset of mitosis inhibits the interaction with p115/USO1, preventing tethering of COPI vesicles and thereby inhibiting transport through the Golgi apparatus during mitosis (By similarity). Also plays a key role in spindle pole assembly and centrosome organization (PubMed:26165940). Promotes the mitotic spindle pole assembly by activating the spindle assembly factor TPX2 to nucleate microtubules around the Golgi and capture them to couple mitotic membranes to the spindle: upon phosphorylation at the onset of mitosis, GOLGA2 interacts with importin-alpha via the nuclear localization signal region, leading to recruit importin-alpha to the Golgi membranes and liberate the spindle assembly factor TPX2 from importin-alpha. TPX2 then activates AURKA kinase and stimulates local microtubule nucleation. Upon filament assembly, nascent microtubules are further captured by GOLGA2, thus linking Golgi membranes to the spindle (PubMed:19242490, PubMed:26165940). Regulates the meiotic spindle pole assembly, probably via the same mechanism (By similarity). Also regulates the centrosome organization (PubMed:18045989, PubMed:19109421). Also required for the Golgi ribbon formation and glycosylation of membrane and secretory proteins (PubMed:16489344, PubMed:17314401)..


Validation Data:

GM130 Rabbit mAb [7lqZ] Images

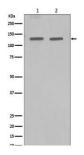

Western blot (SDS PAGE) analysis of extracts from (1) HeLa cell lysate; (2) MCF-7 cell lysate.Using GM130 Rabbit mAb [7lqZ]at dilution of 1:1000 incubated at 4° over night.

View more information on http://naturebios.com

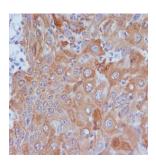

Immunohistochemical analysis of paraffin-embedded human cervix carcinoma, .Using GM130 Rabbit mAb [7IqZ] at dilution of 1:100 incubated at 4 °C over night.Perform heat mediated antigen retrieval before commencing with IHC staining protocol.

IMPORTANT: For western blots, incubate membrane with diluted primary antibody in 1% w/v Milk, 1X TBST at 4°C overnight.

For Research Use Only. Not For Use In Diagnostic Procedures.